




EXOSHAPE® FEMORAL SOFT TISSUE FASTENER


The **EXOSHAPE® FEMORAL SOFT TISSUE FASTENER** is an all-PEEK femoral fixation device that offers aperture soft tissue graft femoral fixation to achieve complete anatomic ACL reconstruction.



Tendon Loop
MAXIMIZES GRAFT LENGTH in tunnel
(At least 20 mm)

COMPRESSED GEOMETRY
designed to reduce the force required for insertion

AN INITIAL SHAPE FOR EASY INSERTION

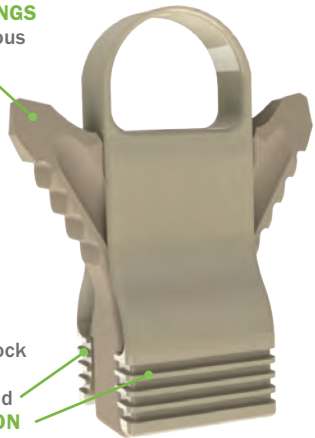

DEPLOYMENT WINGS
lock into cancellous bone

Expandable Block Lock
provides **GRAFT COMPRESSION** and **APERTURE FIXATION**

A FINAL SHAPE FOR SECURE FIXATION

Loaded on a DISPOSABLE DEPLOYMENT GUN



SUTURE CLEAT
to secure loose suture strands

KNOB rotates Implant relative to the Gun. Designed for ergonomics during deployment

DEPLOYMENT TRIGGER
offers easy method to advance Deployment Wings into Tendon Loop

Flexible **CANNULA** for simple navigation inside the joint

Featuring Shape Memory PEEK Altera®



MANUFACTURED WITH
PEEK ALTERA®
SHAPE MEMORY BIOPOLYMER

ACCESSORY INSTRUMENTATION

QuikWhip™

#2 Co-Braid Ultra-High Molecular Weight Polyethylene Surgical Suture

- High Tensile Strength
- Offered in Green and White/Green Co-Braid to help distinguish between graft bundles
- Looped design provides easy method to rapidly stitch graft ends
- Loaded with 51 mm straight round body needle for easy handling



EXOSHAP FEMORAL SOFT TISSUE FASTENER

| Part No. | Description |
|-------------|---|
| 1600-00-008 | ExoShape FEMORAL Cartridge, 8 mm x 21 mm |
| 1600-00-009 | ExoShape FEMORAL Cartridge, 9 mm x 21 mm |
| 1600-00-010 | ExoShape FEMORAL Cartridge, 10 mm x 21 mm |
| 1600-00-011 | ExoShape FEMORAL Cartridge, 11 mm x 21 mm |

SINGLE USE ASSESSORIES

| Part No. | Description |
|-------------|---|
| 1101-03-000 | QuikWhip Green Co-Braid USP #2 Suture, Sterile |
| 1101-04-000 | QuikWhip White/Green Co-Braid USP #2 Suture, Sterile |



CONMED Corporation
525 French Road, Utica, NY 13502

Local 727-392-6464
Toll Free 1-800-237-0169

CONMED.com
customer_service@CONMED.com

For further product information or to arrange a product demonstration, please contact your local CONMED representative or call Customer Service at 800-237-0169. Visit www.CONMED.com.