

Paradigm™ Hip Access Kits*


Featuring innovative products that are designed to provide simpler, safer access to the hip joint with a more efficient, reproducible technique – the Paradigm™ Hip Access Kits provide a convenient, all-in-one solution for hip arthroscopy procedures.




Paradigm™
Tactile Needle*

Paradigm™
Access Needle*


SIGNATURE SERIES™
BLADES AND BURS

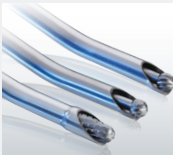

PRESSFT™
SUTURE ANCHORS

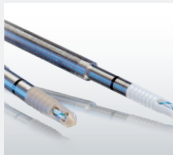


EDGE®
XL PROBES-HIP ARTHROSCOPY

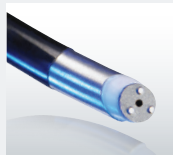

To learn more about these and other innovative products, call **866-4CONMED** or visit **CONMED.com**.

*Not for sale or distribution in all markets.

PARADIGM™
HIP
ACCESS KITS

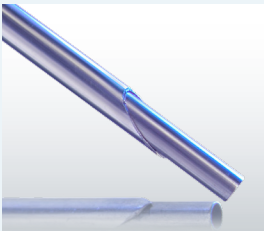

PARADIGM™
TACTILE NEEDLE*

A Paradigm Shift in Accessing the Hip Joint

Introducing an innovative portfolio of products that are designed to be safe, fast, and efficient, allowing surgeons to access the hip joint like never before.

Paradigm™ EZ Switch™ Portal Saver

Say goodbye to the sleds and switching sticks of rigid, motion-restricting cannulas. The innovative EZ Switch™ Portal Saver is flexible to allow a wide range of motion while allowing you to easily and efficiently switch portals.

- **Designed to Save Time** - Easily switch portals without sleds or switching sticks.
- **Flexible Design** - Allows for wide range of motion.
- **Customizable** - Flexible body length (120mm) can be easily trimmed to match patient anatomy.
- **Optimal Access** - Up to 150mm working length is designed to enable easy access to the surgical site.
- **30mm Threads** - Designed to provide strong dermal fixation.

Paradigm™ Tactile Needle*

Engineered for safety, the Tactile Needle is sharp to pierce through skin, muscle, and capsule and then features a spring-loaded blunt obturator for the open joint space that is designed to reduce the risk of damage to the femoral head, cartilage, labrum, and other soft tissues. The needle provides tactile feedback to confirm capsule entry so that surgeons know when they are safely in open space.

Paradigm™ Access Needle*

Designed to more easily penetrate the capsule, Paradigm™ Access Needles feature a sharp 21° bevel as well as a rigid design for optimal maneuvering.



Paradigm™ EZ Switch™
Portal Saver System

Convenient All-In-One Kits

Paradigm™ Hip Access Instruments are available in several kit configurations to meet the needs of your technique and your facility.

Paradigm™ EZ Switch™ Portal Saver System

- (1) Cannulated Obturator
- (1) Reference Tool
- (2) Paradigm™ EZ Switch™ Portal Savers


Paradigm™ Access Kit and EZ Switch™ System*

- (1) Cannulated Obturator
- (1) Reference Tool
- (2) Paradigm™ EZ Switch™ Portal Savers
- (2) Paradigm™ Access Needles*
- (2) Guide Wires

Paradigm™ Tactile Access Kit and EZ Switch™ System*

- (1) Cannulated Obturator
- (1) Reference Tool
- (2) Paradigm™ EZ Switch™ Portal Savers
- (1) Paradigm™ Tactile Needle*
- (1) Paradigm™ Access Needle*
- (2) Guide Wires

PARADIGM™
HIP
ACCESS KITS*

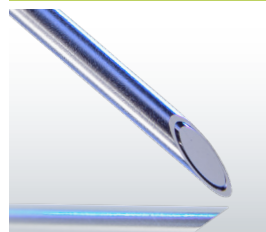

PARADIGM™
ACCESS NEEDLE*

*Experience the Paradigm shift. Take the next step toward simple, safe, and efficient hip arthroscopy. Call **1-866-4CONMED** or visit www.CONMED.com/Paradigm to learn more.*

Ordering Information

	Description	Catalog Number
NEW	Paradigm™ Access	
	Paradigm™ EZ Switch™ Portal Saver System	PEZS01
	Paradigm™ Access Kit and EZ Switch™ System*	PEZS02
	Paradigm™ Tactile Access Kit and EZ Switch™ System*	PEZS03

*Not for sale or distribution in all markets.